





LEGEND LP Tube

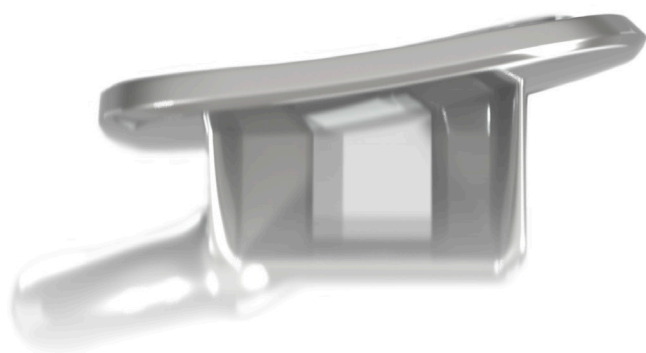


Low Profile Bonding Buccal Tubes
with **Anatomical Micro-Etched Base**

LEGEND LP Tube

LEGEND LP Tube is the latest bonding buccal tube created by TOMY Inc. LEGEND LP Tubes feature a large entry with chamfered edges to facilitate wire insertion along with gingival and occlusal side notches for a secure grip during positioning. LEGEND LP Tubes also include the same successful characteristics found in TOMY's current buccal tubes such as a curved anatomical contour base design, a "V" indent in crown, colour marks for easy quadrant identification and a micro-etched, 80-gauge mesh contour base for superior bonding. Its low profile and lowered mesially positioned hook enhances patient comfort.

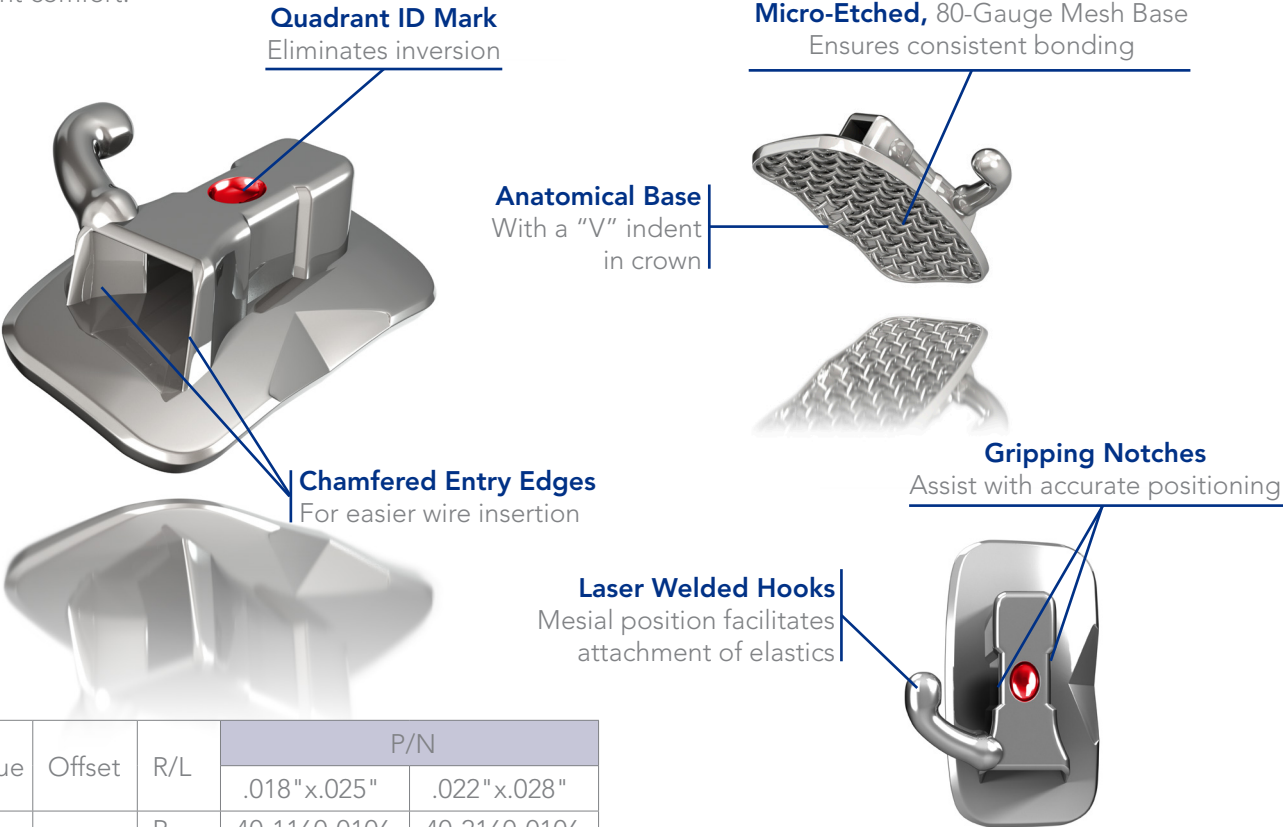

	Torque	Offset	R/L	P/N	
				.018"x.025"	.022"x.028"
U6	-14°	10°	R	40-1160-0106	40-2160-0106
			L	40-1260-0106	40-2260-0106
U6	-14°	0°	R	40-1160-0107	40-2160-0107
			L	40-1260-0107	40-2260-0107
U7	-14°	10°	R	40-1170-0106	40-2170-0106
			L	40-1270-0106	40-2270-0106
U7	-14°	0°	R	40-1170-0107	40-2170-0107
			L	40-1270-0107	40-2270-0107
L6	-20°	0°	R	40-1460-0106	40-2460-0106
			L	40-1360-0106	40-2360-0106
L7	-10°	0°	R	40-1470-0106	40-2470-0106
			L	40-1370-0106	40-2370-0106

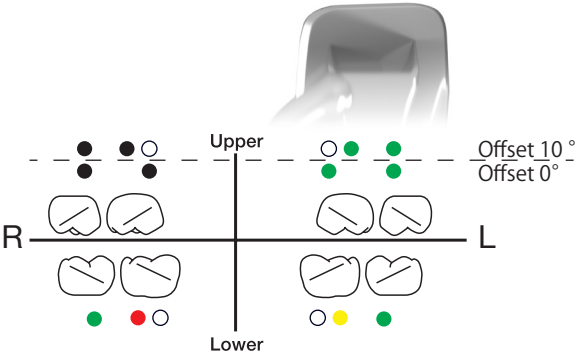

GC Orthodontics Europe GmbH
Harkortstrasse 2
D-58339 Breckerfeld
Germany

