







MEDESY[®]
ITALY



GENERAL CATALOGUE | CATALOGO GENERALE

HIGH QUALITY DENTAL INSTRUMENTS

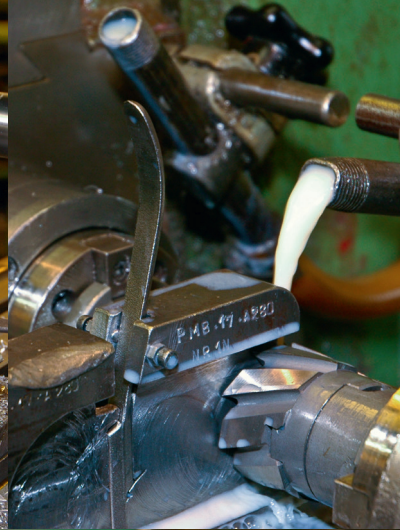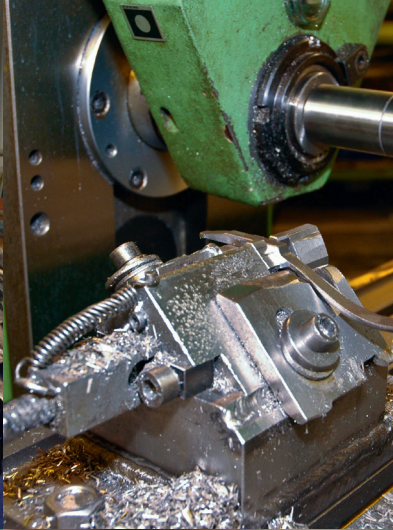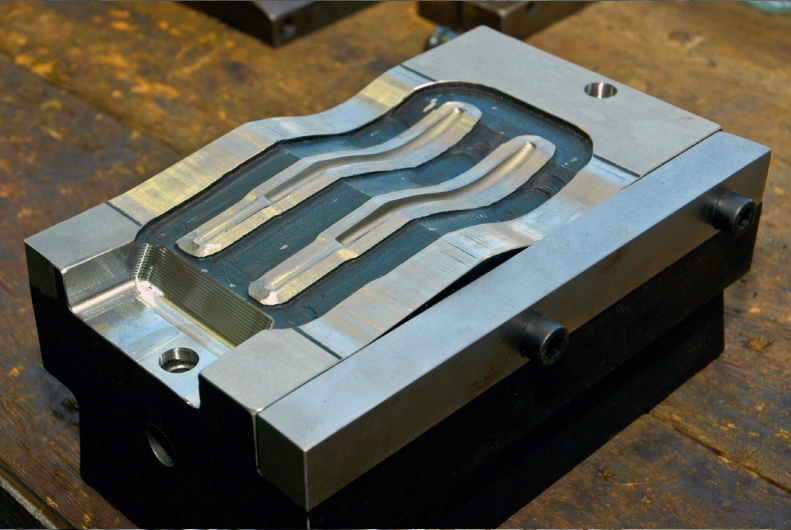

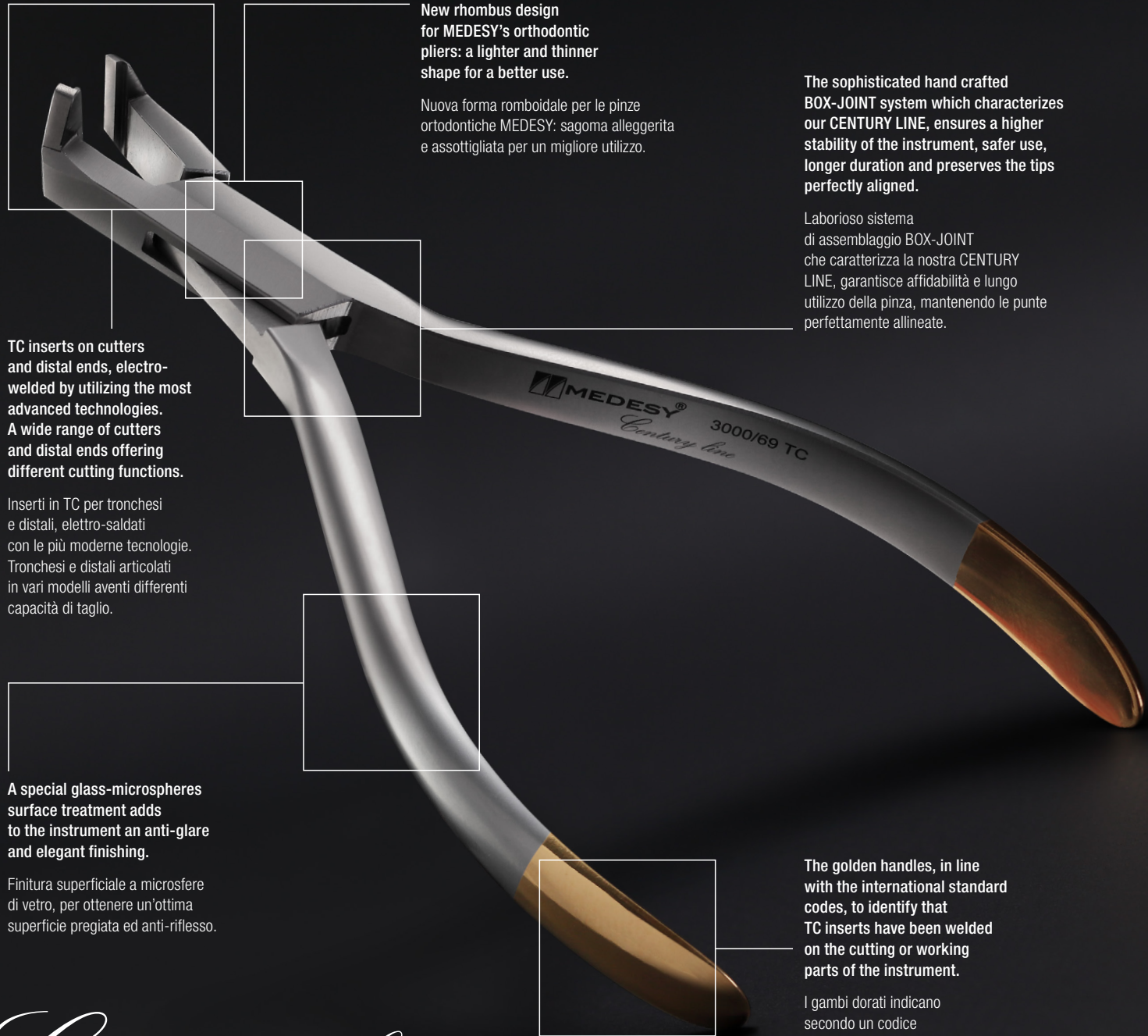

New rhombus design for MEDESYS's orthodontic pliers: a lighter and thinner shape for a better use.

Nuova forma romboidale per le pinze ortodontiche MEDESYS: sagoma alleggerita e assottigliata per un migliore utilizzo.

TC inserts on cutters and distal ends, electro-welded by utilizing the most advanced technologies. A wide range of cutters and distal ends offering different cutting functions.

Inserti in TC per tronchesi e distali, elettro-saldati con le più moderne tecnologie. Tronchesi e distali articolati in vari modelli aventi differenti capacità di taglio.

A special glass-microspheres surface treatment adds to the instrument an anti-glare and elegant finishing.

Finitura superficiale a microsfere di vetro, per ottenere un'ottima superficie pregiata ed anti-riflesso.

The sophisticated hand crafted BOX-JOINT system which characterizes our CENTURY LINE, ensures a higher stability of the instrument, safer use, longer duration and preserves the tips perfectly aligned.

Laborioso sistema di assemblaggio BOX-JOINT che caratterizza la nostra CENTURY LINE, garantisce affidabilità e lungo utilizzo della pinza, mantenendo le punte perfettamente allineate.

The golden handles, in line with the international standard codes, to identify that TC inserts have been welded on the cutting or working parts of the instrument.

I gambi dorati indicano secondo un codice internazionale l'inserimento di placche in TC nella parte lavorante o tagliente.

SPECIFICATIONS ON THE CODES SPECIFICHE SUI CODICI

Code Codice	3000/74 TC
Type Modello	Flush slim
Thickness of wire Spessore del filo	ø max 0,41 × 0,55 mm / .016" × .021"
Length of the instrument Lunghezza dello strumento	(mm 125)

Yellow color for TC instruments
Colore giallo per strumenti in TC

Grey color for standard instruments
Colore grigio per strumenti tradizionali

TC

DISTAL END CUTTERS TRONCHESINI DISTALI CON RITENUTA

3000/69 TC
(1:1 scale / scala 1:1)



TC



3000/69 TC

ø max 0,55 × 0,64 mm / .021" × .025"
(mm 125)

This cutter cuts and holds the portion of wire in excess. Developed for stainless steel and NiTi wires.



TC




3000/69 TC LL

ø max 0,55 × 0,64 mm / .021" × .025"
(mm 150)

Same as 3000/69 TC but with long handle.



TC

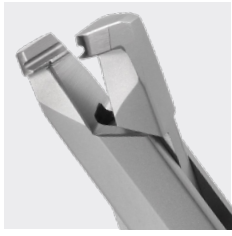


3000/74 TC

Flush slim
ø max 0,41 × 0,55 mm / .016" × .021"
(mm 125)

This cutter cuts and holds wire distally flush to the buccal tube.



TC

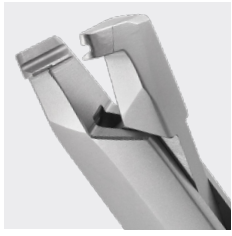

3000/76 TC

Flush
ø max 0,55 × 0,64 mm / .021" × .025"
(mm 125)


Same as 3000/74 TC but with regular head.




WIRE CUTTERS
TRONCHESINI

PIN AND LIGATURE CUTTERS

3000/89 TC
(1:1 scale / scala 1:1)

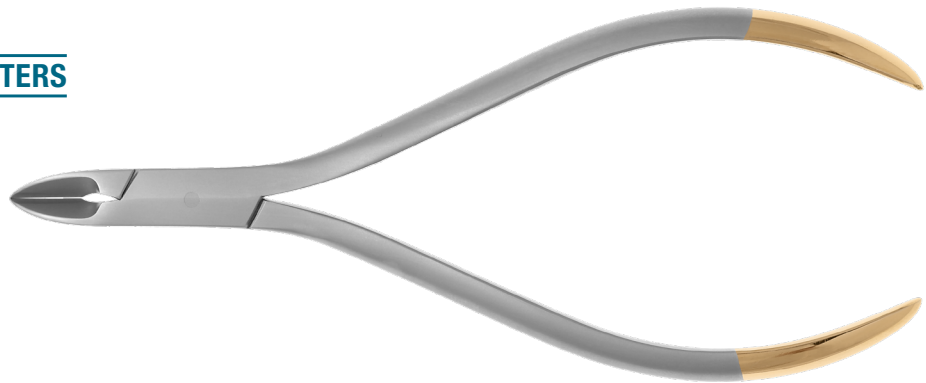

TC

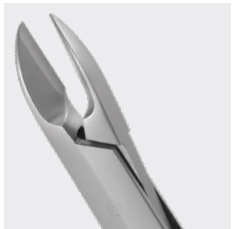


3000/89 TC

ø max 0,3 mm / .012"
(mm 120)

*Delicate cutter, specially designed
for the removal of soft stainless
steel ligatures.*



TC

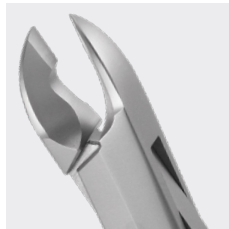


3000/91 TC

ø max 0,3 mm / .012"
(mm 125)

*Delicate cutter with a fine tip for
easier access to difficult areas.
Cut elastics, ligatures, pins and wire.*



TC

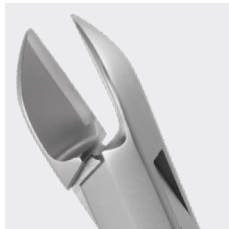


3000/99 TC

ø max 0,4 mm / .016"
(mm 130)

*Medium cutter ensuring accurate
and precision cutting of elastics,
ligatures, pins and wire.*

HARD WIRE CUTTERS

3000/87 TC
(1:1 scale / scala 1:1)

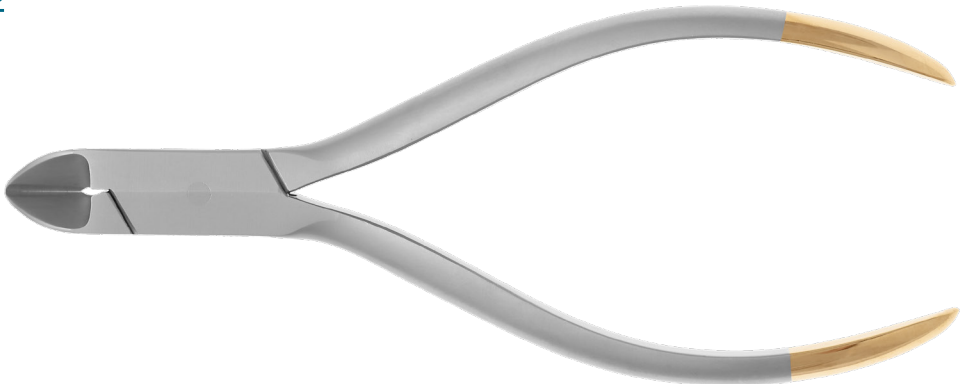

TC



3000/94 TC

ø max 0,4 mm / .016"
(mm 130)

*Medium cutter with rounded
and slim tip. 15° offset tips
help in posterior area.*



TC



3000/87 TC

ø max 0,7 mm / .028"
(mm 130)

*Suitable for cutting archwires
or pins of any size. TC inserts
ensuring longer duration.*



TC



3000/87 TC 15°

ø max 0,7 mm / .028"
(mm 130)

*Same as 87TC, but with a 15°
offset tip enhancing a better
access in difficult areas*



TC



3000/85 TC

ø max 0,3 mm / .012"
(mm 130)

*Standard cutter with slim profile,
allowing an easy access and
better view of the oral cavity.*



TC



3000/85 TC 15°

ø max 0,3 mm / .012"
(mm 130)

Angled 15°



*Same as 85TC, but with a 15°
offset tip enhancing a better
access in difficult areas*

MAX Ø 0.7 mm
.028"

“
MEDESYS IN SPOT
*The orthodontic pliers
are laser engraved with
the maximum cut capacity
(for cutters only) or applicable
wire, in millimeters and inches.*
Le pinze ortodontiche riportano
la capacità massima di taglio
(per i tronchesi) o il diametro massimo
consentito, in millimetri e inches.
”



TC




3000/92 TC

ø max 1,0 mm / .040"
(mm 140)

For front cutting of heavy wires.



TC

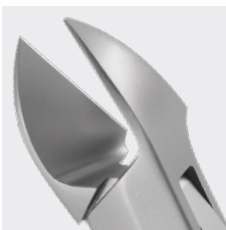


3000/18 TC

ø max 1,0 mm / .040"
(mm 155)

*For an easy cutting of any
size of archwire or pin.
15° offset tips for an easier
access to the posterior area.*

TC

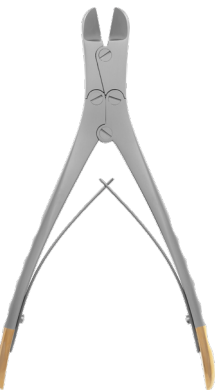

9952/160 TC

ø max 1,5 mm / .060"
(mm 170)



*Double action system for heavy
wires, reduce effort on cutting.*






UTILITY PLIERS
PINZE AUSILIARIE

3000/31
(1:1 scale / scala 1:1)

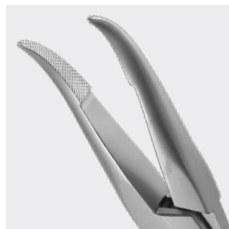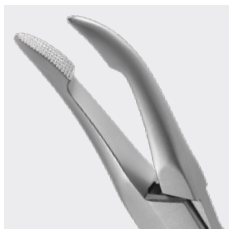

TC

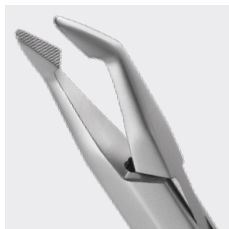

3000/31
Weingart
ø max 0,7 mm / .028"
(mm 145)



Versatile utility plier with very fine serration.

3000/32
Weingart – Slim
ø max 0,5 mm / .020"
(mm 145)

Weingart with fine points.

3000/25 TC
Weingart angled
ø max 0,5 mm / .020"
(mm 140)

Weingart with TC tips, indicated also for placement and removal of archwires, because of the tapered angular design.

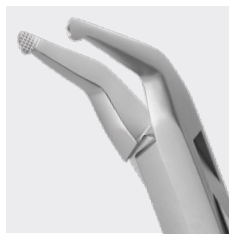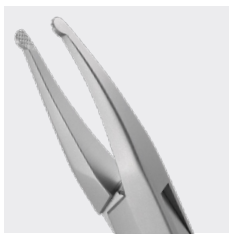

3000/84
How plier – Straight
ø max 0,7 mm / .028"
(mm 140)



Utility plier with round and serrated tips. Ideal for removal and placement of the archwires. Grasping all types of auxiliaries attachments.

3000/72
How plier – Angled
ø max 0,7 mm / .028"
(mm 140)

45° angled tips for convenient access to lingual and others difficult posteriors areas.

WIRE BENDING PLIERS
PINZE PIEGA FILI


3000/39 TC
(1:1 scale / scala 1:1)

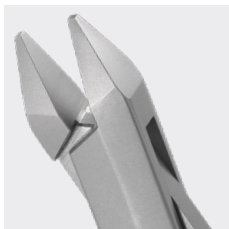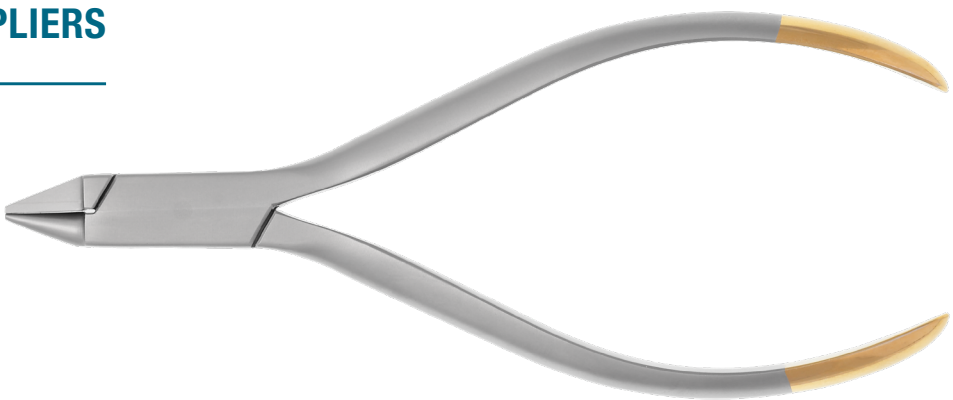

TC



TC

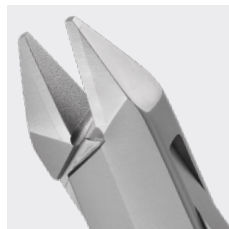

3000/23
Adams
ø max 0,7 mm / .028"
(mm 125)

Squared beaks, clasp and wire bending plier.

3000/23 TC
Adams
ø max 0,7 mm / .028"
(mm 125)


Squared beaks, clasp and wire bending plier, TC inserts.

3000/47 TC
N.64
ø max 1,0 mm / .040"
(mm 140)

Squared beaks, heavy type, TC inserts.



TC

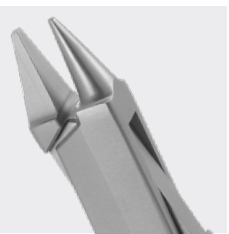

3000/39 TC
Angle N.139
ø max 0,7 mm / .028"
(mm 125)

Bird beak plier, cone and pyramidal tips in TC. Indicated for loops and general bending of round wires.





TC



3000/86 TC
Angle N.138
ø max 1,0 mm / .040"
(mm 130)

Long bird beaks for heavy wire. TC inserts.

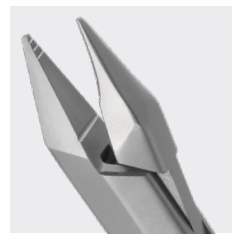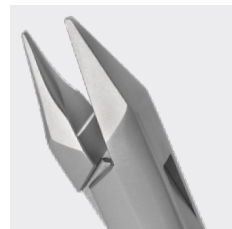

3000/21
ø max 0,7 mm / .028"
(mm 130)

4 grooved bird beaks used to form precise loops in round, squared or rectangular light wire.



3000/54
(1:1 scale / scala 1:1)



TC





3000/40C TC

Tweed
ø max 0,4 mm / .016"
(mm 130)

With cutter

Ideal with round and square wire,
for making consistent uniform loops.
Section tips 1,2-1,8-2,0 mm.

TC

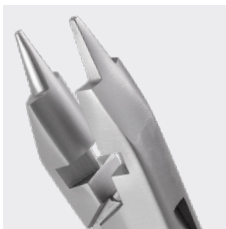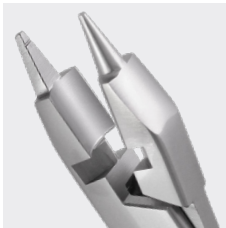

3000/38C TC

ø max 0,6 mm / .024"
(mm 135)

With cutter

Light wire cutting and
contouring plier. For round
wire, to bend small ø loops.
Contouring of loops and arches
up to 0,55 × 0,64 mm.


“ **MEDESY IN SPOT**

Conversion chart
Scala di conversione

mm	Inch
0,3	0,012"
0,4	0,016"
0,5	0,020"
0,6	0,024"
0,7	0,028"
0,8	0,031"
0,9	0,035"
1,0	0,040"


”

3000/51
(1:1 scale / scala 1:1)

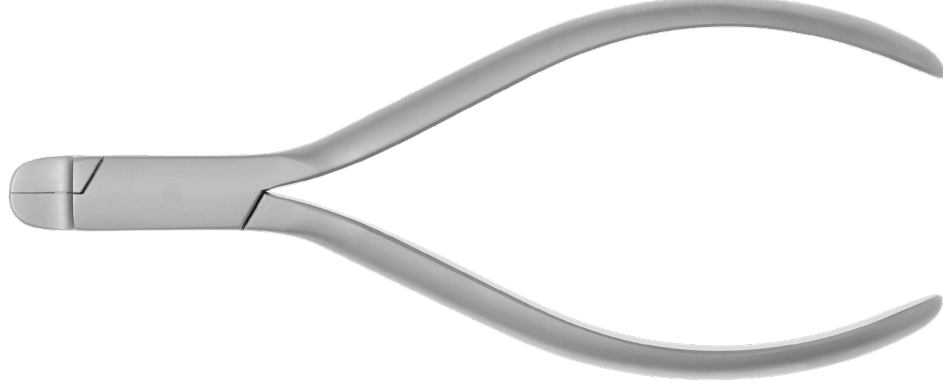

3000/51

Tweed – Short
ø max 0,5 mm / .020"
(mm 125)

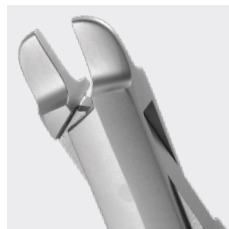

Arch forming plier – Short.
Ideal for torquing and bending
squared or rectangular wire.
Parallel blades prevents
crushing of wires.



3000/52

Tweed – Long
ø max 0,7 mm / .028"
(mm 130)



Arch forming plier – Long.
As 3000/51 but for bigger wires size.



3000/30

Schwarz
ø max 0,7 mm / .028"
(mm 130)




Arrow clasp bending plier.



3000/57

Aderer – Short
ø max 0,6 mm / .024"
(mm 120)

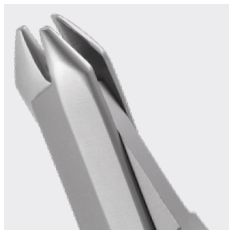


3 jaw plier – Short.
Slightly rounded tips,
for contouring and bending
without damaging the archwire.



3000/54

Aderer – Long
ø max 0,7 mm / .028"
(mm 130)

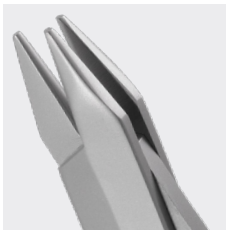


3 jaw plier – Long.
Slightly rounded tips,
for contouring and bending
without damaging the archwire.



3000/20

Jarabak
ø max 0,7 mm / .028"
(mm 140)

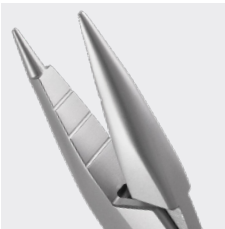

Light wire bending plier.
Versatile in different uses,
as bending or closing loops.
One beak partially serrated
for safe grip.



3000/34

De La Rosa
ø max 0,5 mm / .020"
(mm 125)




Contours and forms archwires,
of light wire without torquing.



3000/35

De La Rosa
ø max 0,5 mm / .020"
(mm 125)

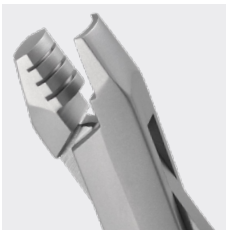

Same as 3000/34 with additional
grooves measuring respectively
0,41-0,46-0,56 mm.





3000/58

ø max 0,7 mm / .028"
(mm 125)



Lingual arch forming plier for
fast forming all lingual retainers.

3000/98
(1:1 scale / scala 1:1)

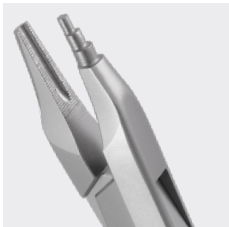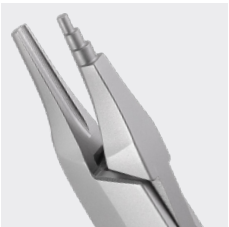

3000/98
Tweed – Short
ø max 0,5 mm / .020"
(mm 135)


Designed to form precise omega loops. Long grooved beak internally smoothed. Short type.

3000/95
Tweed – Short
ø max 0,5 mm / .020"
(mm 135)

Designed to form precise omega loops. Short grooved beak internally serrated. Short type.

3000/55
Young
ø max 0,7 mm / .028"
(mm 130)

Ideal for labial arches and face bow. Universal plier for different size loops.

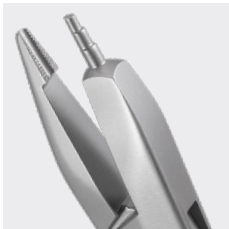

3000/96
Tweed – Long
ø max 0,5 mm / .020"
(mm 145)

Designed to form precise omega loops. Long grooved beak internally smoothed. Long type.

3000/97
Tweed – O'Brein
ø max 0,7 mm / .028"
(mm 140)

Graduated cone to form precise omega and closing loops. Micro serrated beaks to prevent wire slippage.

3000/81
(1:1 scale / scala 1:1)



• ø 0,7 mm



3000/53
Optical
ø max 0,7 mm / .028"
(mm 140)

Concave/convex loop forming plier, small-medium.

• ø 1,2 mm



3000/79
Optical
ø max 0,7 mm / .028"
(mm 135)

Concave/convex loop forming plier, medium-large.



3000/70
Hollow chop
ø max 0,7 mm / .028"
(mm 140)

Contours and form arches of all sizes.

•• ø 0,5 mm



3000/81
Long and delicate beaks plier.
ø max 0,5 mm / .020"
(mm 135)

••



3000/28
2 round beaks long.
ø max 0,5 mm / .020"
(mm 140)

MEDESy IN SPOT

The new lubricant oil for orthodontic instruments. If used regularly it ensures a long product lifespan and optimal performances: facilitates the smooth and perfect function for all hinged instruments and helps to prevent freezing from corrosion and sediments.


Nuovo olio lubrificante per pinze ortodontiche. Se usato regolarmente prolunga la durata e l'efficienza dello strumento: l'olio favorisce il movimento scorrevole e delicato degli strumenti a snodo, proteggendo al tempo stesso da corrosione e sedimenti.

4160
12 ml
Tip length 40 mm
Lunghezza punta 40 mm

4155
12 ml
Tip length 7 mm
Lunghezza punta 7 mm



3000/63
(1:1 scale / scala 1:1)

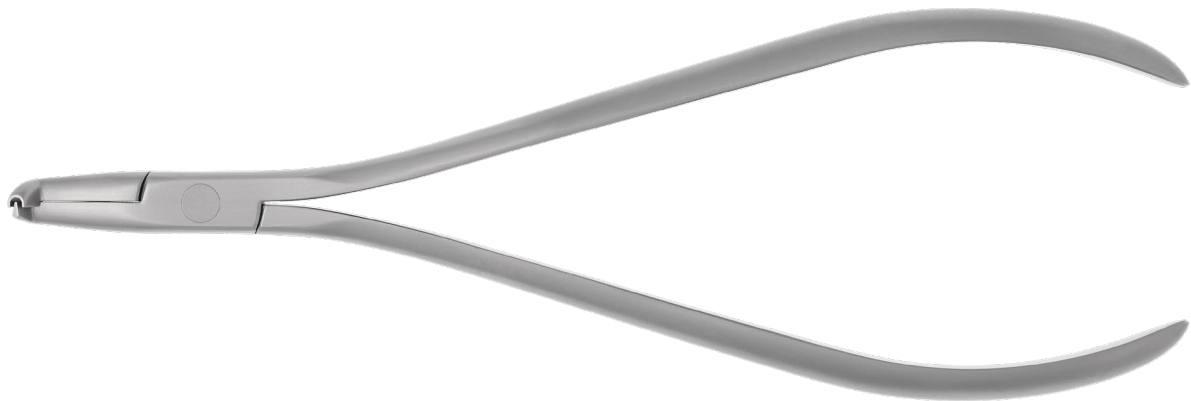

3000/61
(1:1 scale / scala 1:1)



3000/80
Nance
ø max 0,7 mm / .028"
(mm 125)

Loop bending plier.
Four step tips (3, 4, 5, 6 mm)
for loops forming.
Round and rectangular wire.

3000/63
Cinch back
ø max 0,55 × 0,64 mm
ø max .021" × .025"
(mm 155)

Cinch back utility plier.
For closing arches or cinch down a wire.
Use on round or rectangular wire.



3000/59
Johnson
ø max 0,7 mm / .028"
(mm 135)

Crown and band contouring plier.

3000/61
Reynolds
ø max 0,7 mm / .028"
(mm 130)

Crown and band contouring plier.

3000/29
Schwarz
ø max 0,7 mm / .028"
(mm 125)

Arrow clasp bending plier.



3000/43
ø max 0,5 mm / .020"
(mm 140)

Closing loop plier.
Different loops in only one plier.

3000/60
ø max 1,4 mm / .055"
(mm 130)

3 prong plier. For hard wires.
Round tips prevents any
damage of archwire.

3000/56
Nance
ø max 0,7 mm / .028"
(mm 135)

Loop closing plier.
For bending of stop loops and
ideal for Adams clasps.



3000/15
Step plier
ø max 0,55 × 0,64 mm
ø max .021" × .025"
(mm 140)

Double side head to create
steps in archwire.

3000/171
V-Stop
ø max 0,55 × 0,64 mm
ø max .021" × .025"
(mm 130)

For bending archwire in V-form,
in order to avoid any movement.





3000/49
(1:1 scale / scala 1:1)



3000/83
Bracket remover – Straight
(mm 135)



3000/82
Bracket remover – Angled
(mm 130)

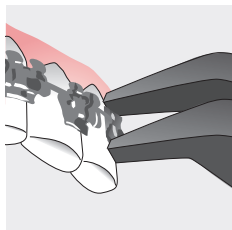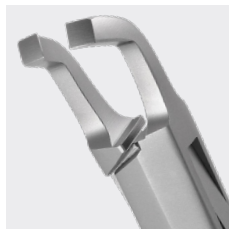


Pliers to remove the brackets: the extremely sharp tips allow an easy access under the flange of the bracket and facilitate its detachment.



Pinze per rimozione attacchi: le punte affilate permettono di entrare sotto la flangia dell'attacco e ne consentono il distacco con estrema facilità.



3000/49
Posterior band remover – Short
(mm 135)

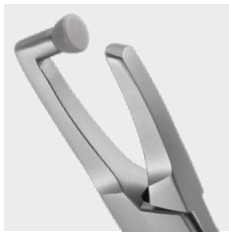

3000/50
Posterior band remover – Long
(mm 135)

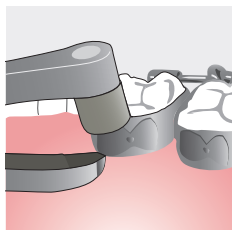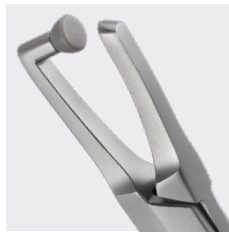

Pliers to remove the bands available in two sizes, short and long. The rubber insert fully autoclavable sticks perfectly to the crown and allows an optimal position of the tip during the removal of the band.

Pinze per la rimozione di bande disponibili in 2 misure, lunga e corta. Il gommino sterilizzabile, rimane perfettamente in contatto con la corona, permettendo una ottimale posizione della punta, per la rimozione della banda.




3000/PR5
Tips spare parts, screwable
Gommini di ricambio, avvitabili

3000/45
(1:1 scale / scala 1:1)



3000/45
Adhesive remover
(mm 140)

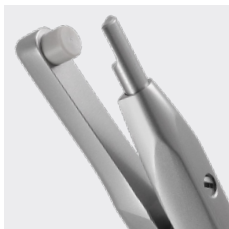

Used to remove composite material after debonding attachments.



3000/PR8
Tips spare parts, screwable
Gommini di ricambio, avvitabili




3000/62
Crimpable archwire hook
(mm 130)



For ball hooks and auxiliaries.



3000/36
Crimping plier
(mm 125)

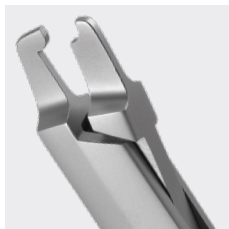

For crimping bands or hooks, for placing stops and tubing on archwires.





LIGATURE FORCEPS
PINZE PER LEGATURE




1851
(1:1 scale / scala 1:1)

Smooth handle
for Smaha forceps
Manico liscio
per le pinze Smaha




1851 (mm 140)
Smaha


1852 (mm 140)
Smaha
Curved
Curvo

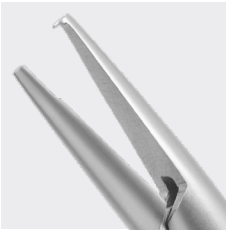

2825 (mm 140)
Plier to place modules
Pinza per moduli

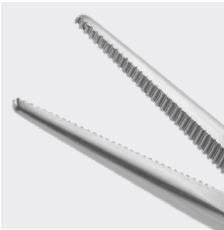


1519 (mm 125)
Halstead Mosquito

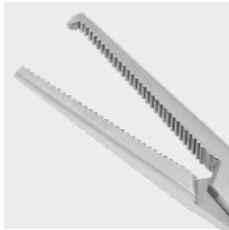

2824 (mm 120)
Forceps for placing dentalastics
Pinza per applicare i dentalastics

2825/ORTO
(1:1 scale / scala 1:1)

Grooved handle
for Mathieus
Manico zigrinato
per Mathieu

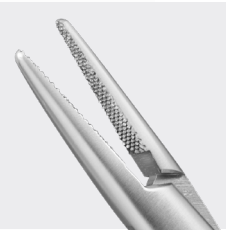


1850 (mm 140)
Mini Mathieu




2825/ORTO (mm 140)
Mini Mathieu smoothed
Smoothed
Liscio



2826 (mm 165)
Ligature plier with locker
Pinza per legature con fermo

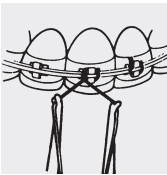

2822 (mm 150)
Steiner
Ligature plier, for applying and stretching
stainless steel ligatures
Pinza per legature, per applicare e tendere le legature

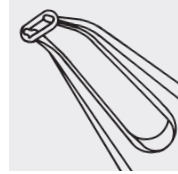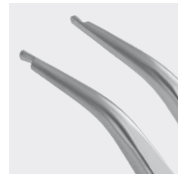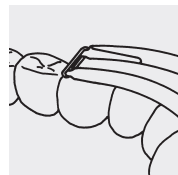


2823 (mm 150)
Separating plier, for an easy placement
of separating elastics
Pinza per separazione, per inserire facilmente gli elastici



ORTHODONTICS / ORTODONZIA

(1:1 scale / scala 1:1)

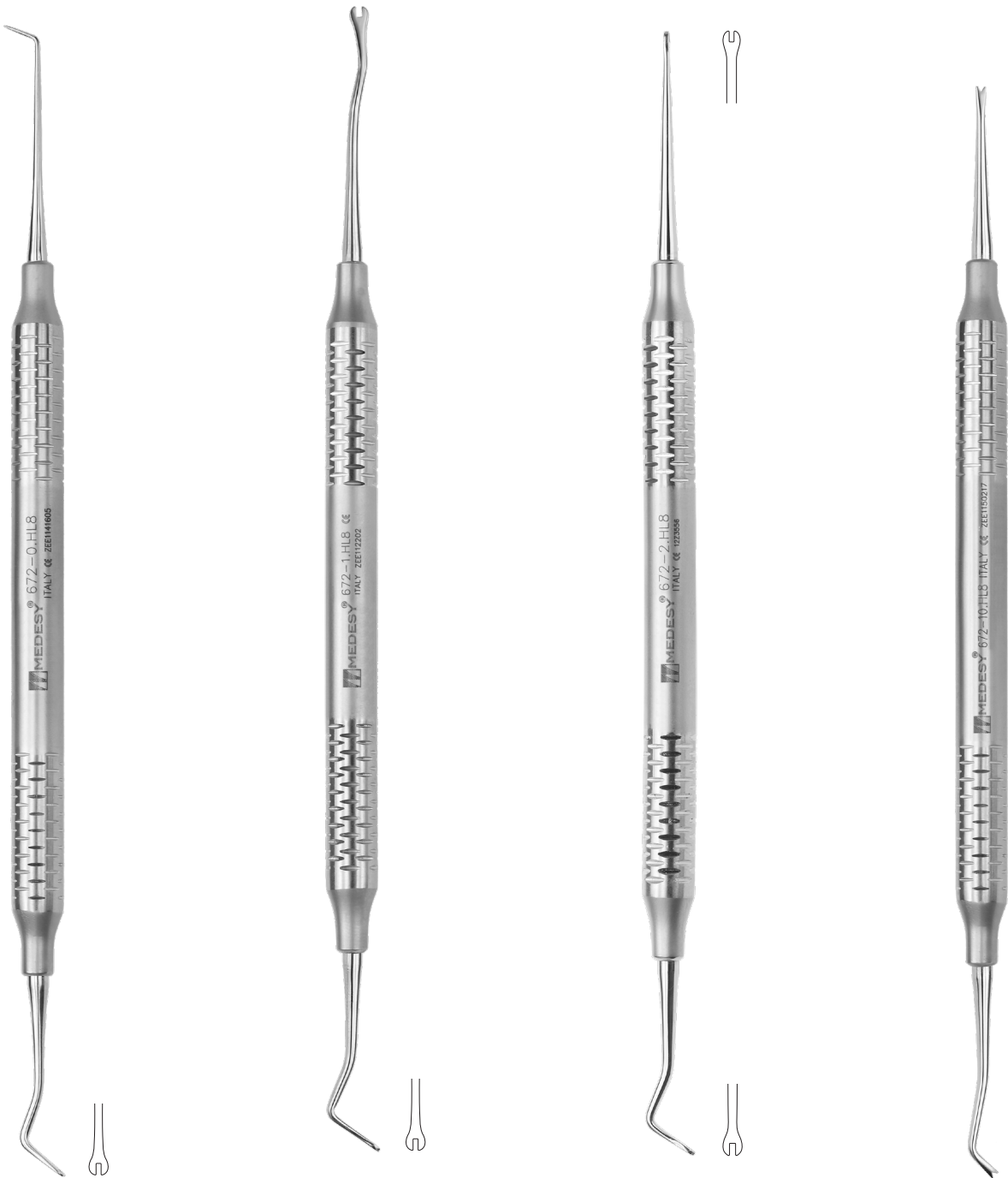

672/0 (ø mm 6)
672/0.HL8 (ø mm 8)
672/0.HL10 (ø mm 10)

Ligature director N. 0
Direzionatore per legature N. 0

672/1 (ø mm 6)
672/1.HL8 (ø mm 8)
672/1.HL10 (ø mm 10)

Ligature director
Direzionatore per legature


672/2 (ø mm 6)
672/2.HL8 (ø mm 8)
672/2.HL10 (ø mm 10)

Ligature director
Direzionatore per legature

672/10 (ø mm 6)
672/10.HL8 (ø mm 8)
672/10.HL10 (ø mm 10)

Ligature director MINI
Direzionatore per legature MINI

(1:1 scale / scala 1:1)

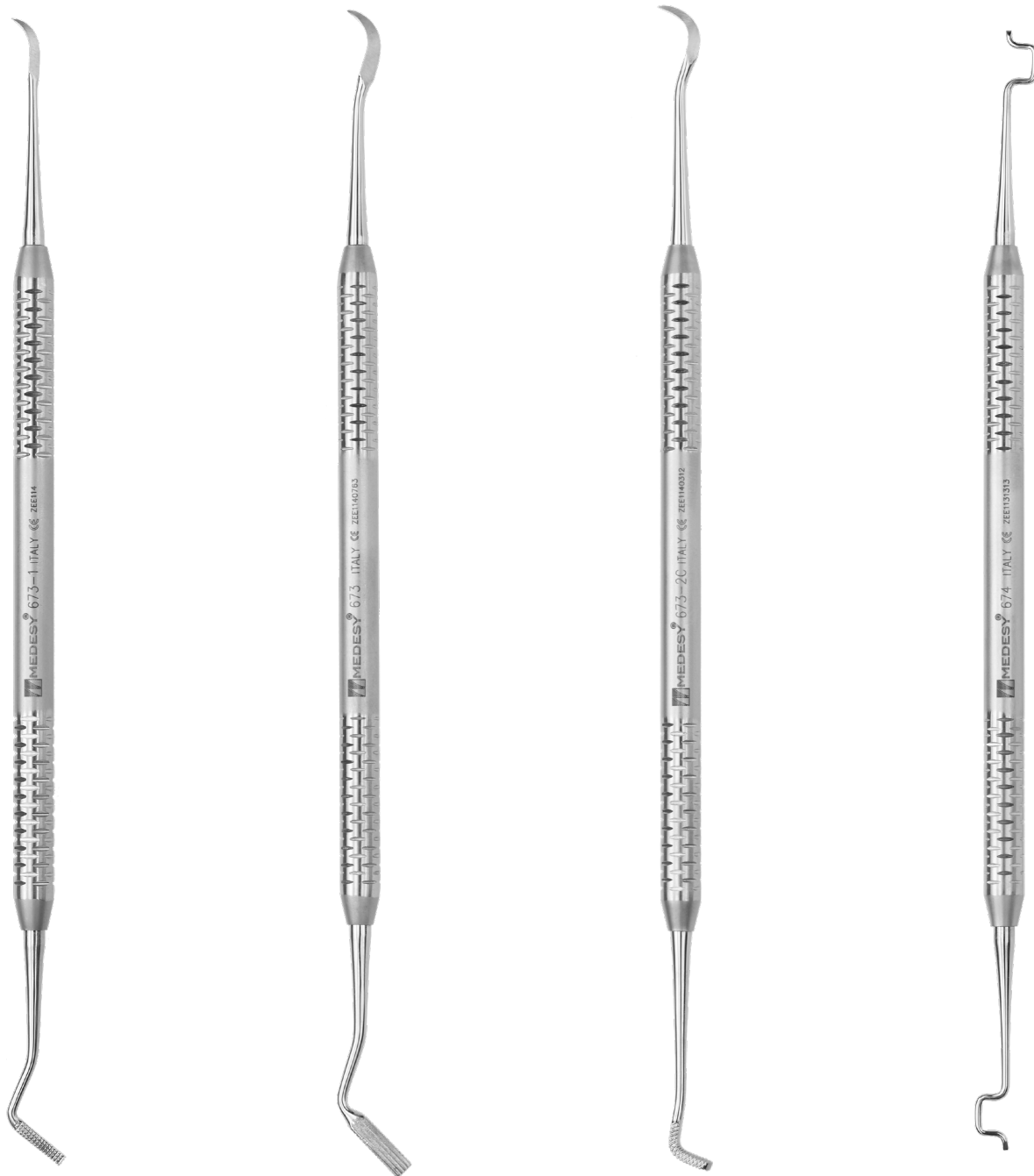

673/1 (ø mm 6)
673/1.HL8 (ø mm 8)
673/1.HL10 (ø mm 10)

Band pusher / scaler MINI
Spingibande / scaler MINI

673 (ø mm 6)
673.HL8 (ø mm 8)
673.HL10 (ø mm 10)


Band pusher / scaler
Spingibande / scaler

673/2C (ø mm 6)
673/2C.HL8 (ø mm 8)
673/2C.HL10 (ø mm 10)

Band pusher / scaler
Spingibande / scaler

674 (ø mm 6)
674.HL8 (ø mm 8)
674.HL10 (ø mm 10)

Elastic ligature applicator
Applicatore per elastici

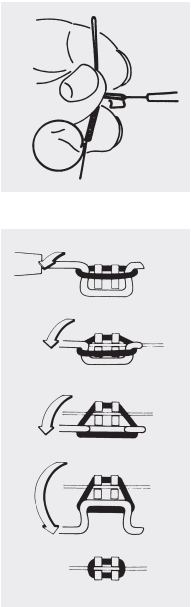

DIRECT BOND BRACKET TWEEZER
PINZETTE PER L'ESATTO POSIZIONAMENTO DEGLI ATTACCHI



1006 (mm 135)
With blade
Con lama


1008 (mm 140)

1009 (mm 140)

1130 (mm 140)

1132 (mm 120)

1134 (mm 140)
Buccal tube tweezer
Pinzetta ausiliaria buccale

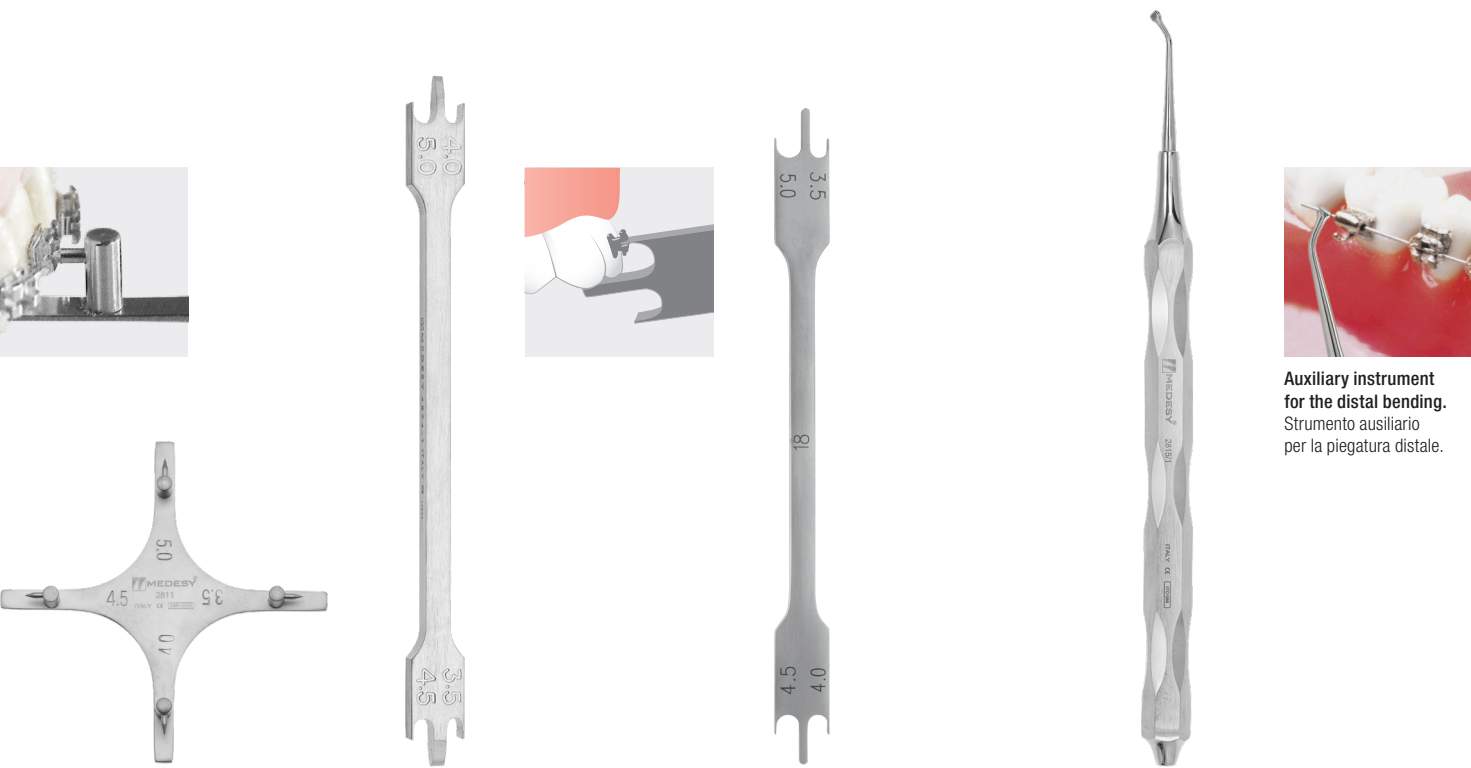

2811
Positioning gauge
Posizionatore per brackets
(mm 65×65)

4874/1 (per slot 018)
4874/2 (per slot 022)
Bracket height Gauge – Aluminium
Altimetro per bracket – Alluminio
(mm 150)

4876/1 (per slot 018)
4876/2 (per slot 022)
Bracket height Gauge – Stainless steel
Altimetro per bracket – Acciaio
(mm 130)

2815/1 (Hole dimension: $\varnothing$ mm 0,75 – 030")
2815/2 (Hole dimension: $\varnothing$ mm 0,55 – 022")
(mm 160)



Auxiliary instrument
for the distal bending.
Strumento ausiliario
per la piegatura distale.



4870 (mm 150)
Mershon
Band pusher
Spingibande

4873 (mm 140)
Molar band seater, adjustable
Posizionatore per bande, regolabile

4871 (mm 145)
Band seater, autoclavable teflon handle
(available in different colours)
Posizionatore per bande, manico in teflon
autoclavabile (disponibile in diversi colori)

2817 (mm 165)
Positioner, aluminium
Posizionatore, alluminio

CROWN AND GOLD SCISSORS



3550 (mm 110)
Beebee
Saw edge, straight
Dentata, retta

3551 (mm 110)
Beebee
Saw edge, curved
Dentata, curva

WIRE SCISSORS



3522 (mm 125)
Universal, saw edge
Universale, dentata

3522/TC (mm 125)
Universal
Universale



ORTHODONTIC KITS
KIT ORTODONTICI

1680/1

ORTHO CLASSIC

Content / Contenuto:

- | | |
|------------|------------------------------|
| 1850 | Mini Mathieu |
| 3000/31 | Weingart plier |
| 3000/69 TC | Distal end n. 69 |
| 3000/85 TC | Cutter n. 85 |
| 4905/8 | Handle with mirror #4 |
| 672/1.HL8 | Ligature director Tucker n.1 |
| 673.HL8 | Band pusher / scaler |
| 981/19 | Gammatrix tray Quatri |



“ MEDESIV IN SPOT

**Content can be changed
under request**

*Il contenuto può essere
cambiato a richiesta*

”

1680/2

ORTHO ADVANCED

Content / Contenuto:

- | | |
|------------|------------------------------|
| 1006 | Bracket tweezer |
| 1850 | Mini Mathieu |
| 3000/23 TC | Adams plier |
| 3000/31 | Weingart plier |
| 3000/69 TC | Distal end n. 69 |
| 3000/85 TC | Cutter n.85 |
| 4905/8 | Handle with mirror #4 |
| 672/1.HL8 | Ligature director Tucker n.1 |
| 673.HL8 | Band pusher / scaler |
| 674.HL8 | Elastic ligature applicator |
| 981/19 | Gammatrix tray Quatri |

1680/3

ORTHO BAND

Content / Contenuto:

- | | |
|------------|-----------------------|
| 3000/31 | Weingart plier |
| 3000/49 | Band remover short |
| 3000/50 | Band remover long |
| 4870 | Band pusher |
| 4873 | Molar band seater |
| 4905/8 | Handle with mirror #4 |
| 673.HL8 | Band pusher / scaler |
| 673/2C.HL8 | Band pusher / scaler |
| 981/19 | Gammatrix tray Quatri |



1680/4

ORTHO BRACKET

Content / Contenuto:


- | | |
|---------|-----------------------|
| 1006 | Bracket tweezer |
| 1008 | Bracket tweezer |
| 1009 | Bracket tweezer |
| 2817 | Positioner |
| 3000/31 | Weingart plier |
| 3000/82 | Bracket remover |
| 3000/83 | Bracket remover |
| 4876/1 | Bracket height gauge |
| 4876/2 | Bracket height gauge |
| 4905/8 | Handle with mirror #4 |
| 981/19 | Gammatrix tray Quatri |




ORTHODONTICS / ORTODONZIA




6141
240 × 120 mm
Large rack
Porta pinze lungo

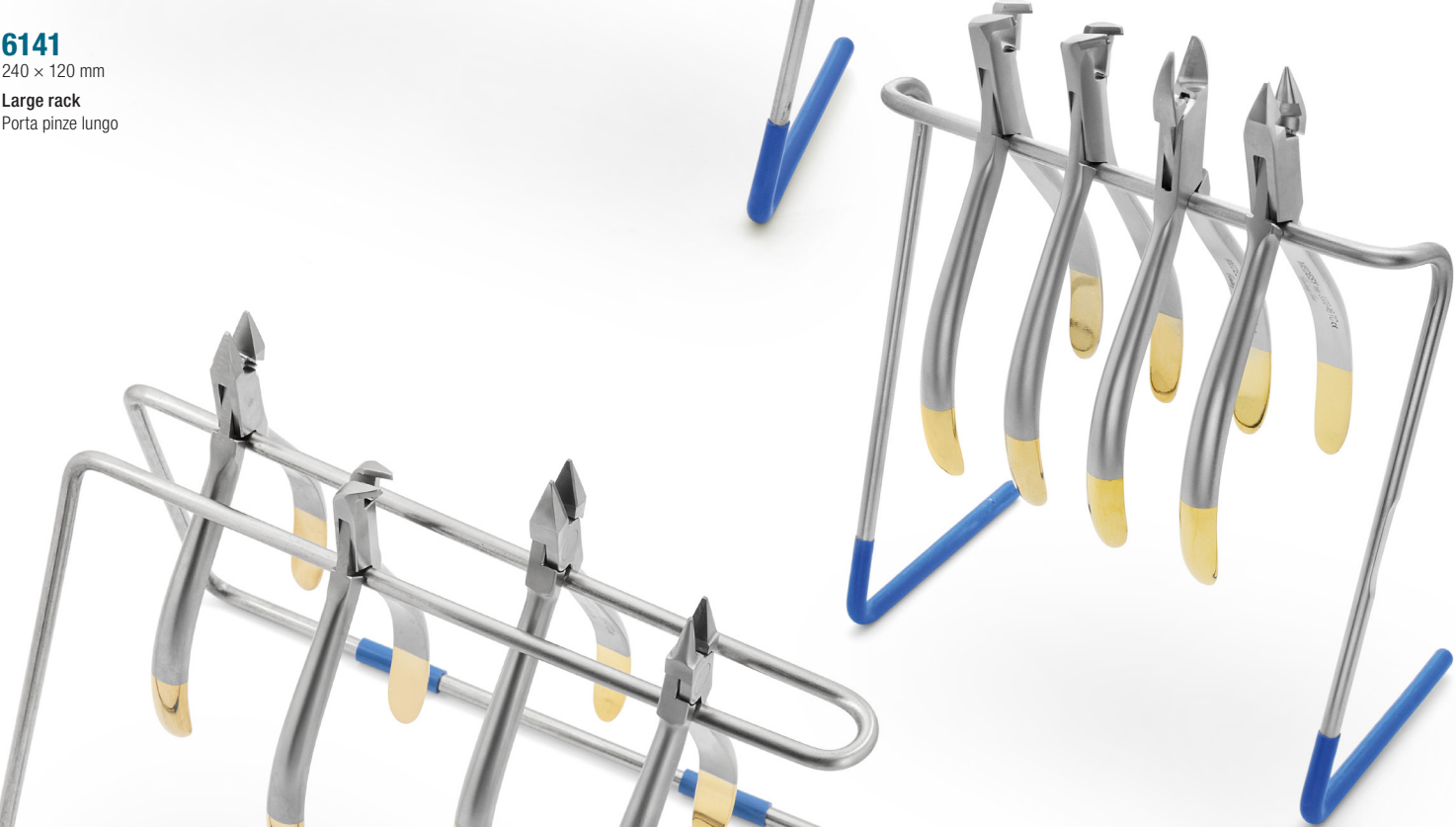

6140
150 × 120 mm
Small rack
Porta pinze corto

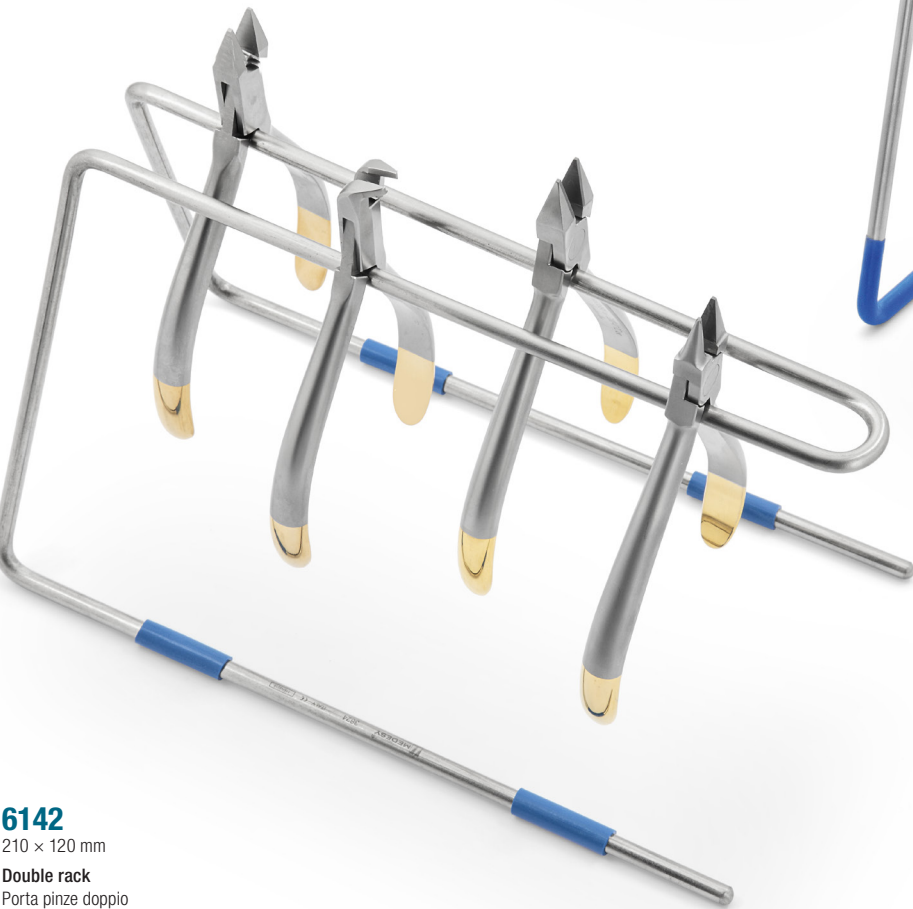

6142
210 × 120 mm
Double rack
Porta pinze doppio

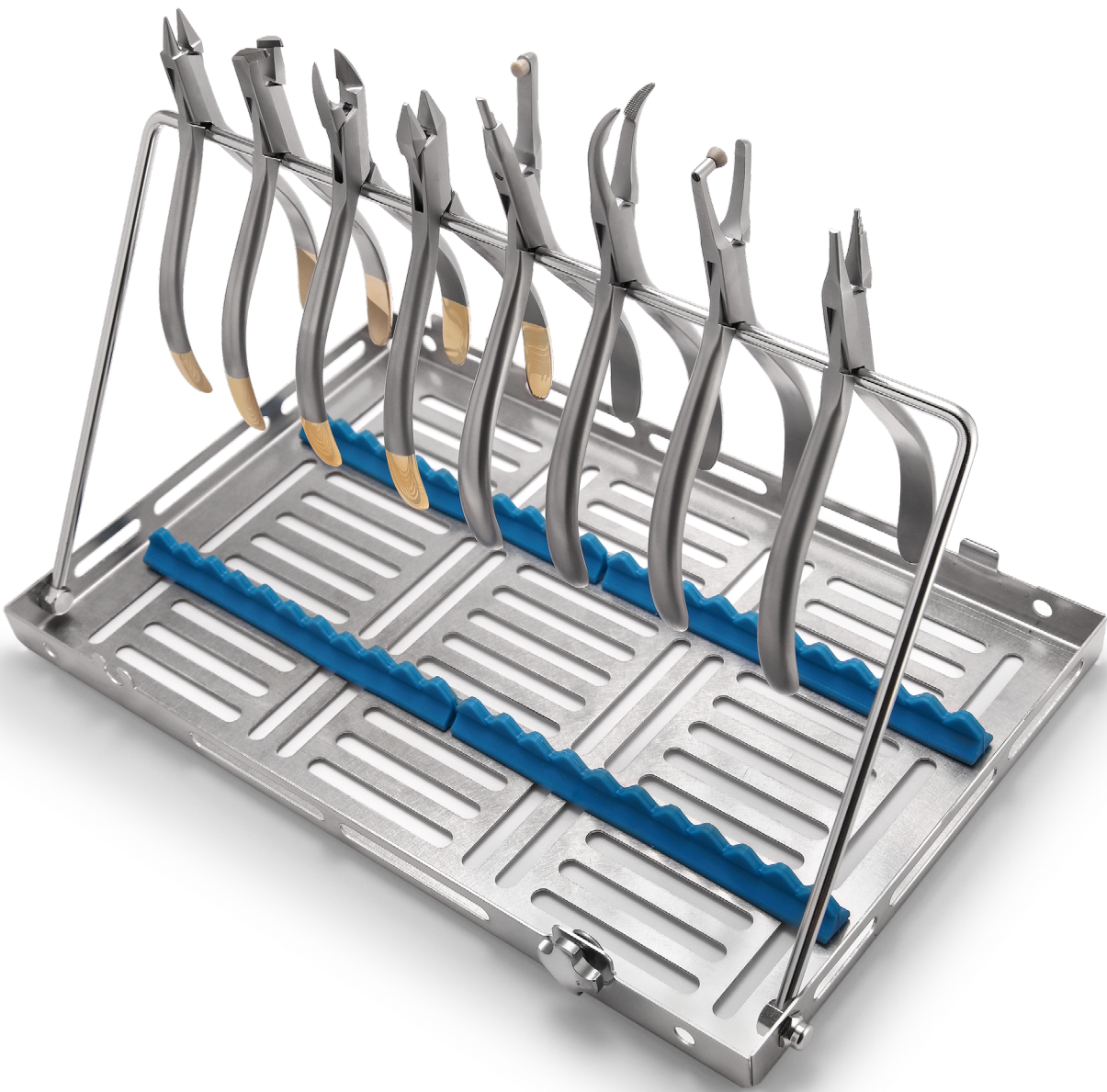


MEDESY offers a very practical solution for handling your pliers during cleaning and sterilization sections. Our GAMMAFIX trays helps on organizing your pliers in practical set ups and handle them safely. Just some of the specific benefits of our GAMMAFIX trays: inserts made of special non-toxic silicon and open corners to facilitate the outflow of residuals and water and therefore reduce the drying time. More info on www.medesy.it

MEDESY presenta un'ottimale soluzione per maneggiare le pinze ortodontiche durante il processo di pulizia e sterilizzazione. I vassoi della serie GAMMAFIX risultano estremamente pratici e versatili per contenere le pinze ortodontiche in tutta sicurezza. I vassoi GAMMAFIX si distinguono per la loro praticità d'uso ed importanti accorgimenti, tra i quali il silicone atossico e angoli aperti atti a defluire acqua e residui. Maggiori informazioni su www.medesy.it



981/ORTO
Based on 981/20
Gammaxif (Cinc)
Basato su 981/20
Gammaxif (Cinc)

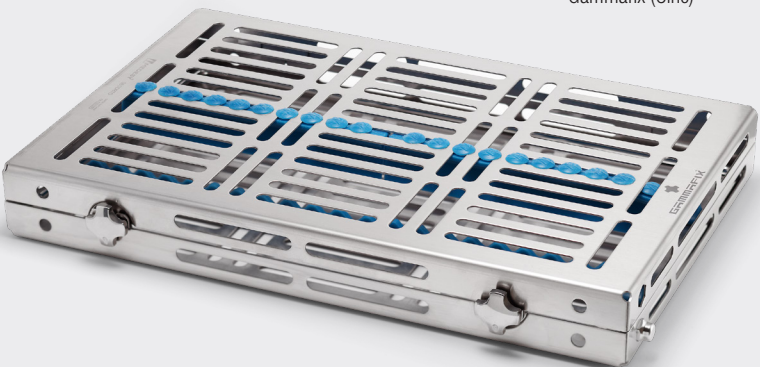



MEDESY presents the practical and elegant display case solution for orthodontic pliers.
Select the pliers according to the needs and carry them, very easy and comfortable.
The perfect choice for congresses and exhibitions!

MEDESY presenta la pratica ed elegante valigetta per pinze ortodontiche.
Seleziona le pinze in funzione alla necessità per trasportarle in modo comodo e semplice.
Ideali per congressi e fiere!



2918/0
cm 55 x 22 x 8
Exposer
Espositore

“

MEDESY IN SPOT

*Content can be changed
under request*

*Il contenuto può essere
cambiato a richiesta*

”

**MEDESY SRL**

Viale dell'Industria, 1 – Industrial area
33085 Maniago, PN – Italy
Tel. (+39) 0427.72786
Fax (+39) 0427.71541
www.medesy.it – info@medesy.it